



Groove™
SERIES

The Power ³

Just imagine...

OPTIONS AND ACCESSORIES





The Power³

OPTIONS AND ACCESSORIES

Just Imagine...

A power mobility solution that is truly tailored around the needs of the client!

A choice of seating tailored to specific clinical needs, a choice of Power bases to suit the environment and an intelligent drive system which is customised to both the user and their surrounding environment.

Meet the Groove series from Quickie, the first series of Power Chairs that breaks the mould with design, modularity, performance and the ability to meet any user's needs without compromise.

This brochure outlines the features and options that the Groove series offers. You will also find a wide range of accessories that complement the Quickie Groove, all of which can be ordered and installed on chair construction or added to the existing chair.

Introduction

Please note options displayed in this catalogue may not be available in your country. All options are subject to compatibility with other options you may have chosen. Always refer to the Quickie Groove Order Form for option availability and compatibility.

Please note the option codes in this brochure are NOT part numbers and cannot be ordered. Option codes relate to the order form and are for reference purposes only.

Please refer to the Parts Manual for Part Numbers.

OPTIONS AND ACCESSORIES

CONTENTS

SECTION 1 - SEATING

Rehab Seating	4
Jay Seating	6
Comfort Seating	8
Perfect Fit Seating	12
Recaro Seating	14

SECTION 2 - POSITIONING

Backrest Recline Options - Standard Rehab & Comfort Seating	16
Backrest Recline Options - Perfect Fit	17
Seat to Base Interfaces	18

SECTION 3 - ARMRESTS & SKIRT GUARDS

Armrests	20
Sideguards	21
Armpads & Handpads	22
Comfort Seating Armpads	23

SECTION 4 - LEGRESTS

Swing away legrests	24
Footplates	27
Legrest Accessories	28

SECTION 5 - Quickie Electronics CONTROL SYSTEMS

QC Controls	30
QR Controls	31
Direct Actuator Control Switches	32
Control Mountings	34
Control Boxes	38
Special Controls	41
Head Control Mountings	42

SECTION 6 - WHEELS & CASTORS 44

SECTION 7 - BATTERIES & CHARGERS 46

SECTION 8 - ACCESSORIES 48

SEATING

With versatile and adaptable seating options, the range of Quickie Groove seating accessories is extensive. The comprehensive selection of cushions, backs and positioning aids means the Quickie Groove is ready to meet all client demands whatever their lifestyle or clinical requirements.

Upgrading, retro fitting or reconfiguring has never been easier than with the Quickie Groove, so as a client's needs change so does their Groove.

REHAB SEATING

REHo20100

Solid seat pan

The solid seat pan is the base for the seat cushion and is available in the following widths to suit all seat depths:

36 cm (14")
41 cm (16")
46 cm (18")
51 cm (20")
56 cm (22")



REHo20200

5 cm (2") standard foam cushion

Our standard foam cushion is available in sizes to suit all seat widths and depths



REHo20101

Contoured seat cushion

The Standard Rehab contoured seat cushion provides clients with the added support around the knees and legs to aid positioning.

Available in sizes to suit all widths and depths



OPTIONS AND ACCESSORIES



Standard back sling

Our tried and tested Quickie back sling is the standard feature on Standard Rehab seating, padded for extra comfort and ideal for clients with low clinical requirements.

Available in sizes to suit backrest widths 36cm - 56cm (14" - 18") & 38cm - 51cm (15" - 20") Tall

REHo30000



Solid backrest

Padded for extra comfort the detachable solid backrest is ideally suited for clients requiring extra support and stability.

Available in the following sizes:

36cm (14") wide - 38cm - 51cm (15 - 17") Tall

41cm - 56cm (16" - 22") wide - 51cm (20") Tall

REHo30001



Optima backrest

The optima backrest provides excellent support and flexibility through tension adjustable straps which can be adjusted and shaped to the client's back for greater comfort and better posture.

Available in sizes to suit backrest widths 36cm - 46cm (14" - 18").

REHo30002

JAY® SEATING SYSTEM

For clients with varying posture, skin and clinical requirements, Jay seats and backrests provide the support and comfort not provided by other seating products.

Jay seats stabilise, support and position the body while the range of Jay backs provide posterior lateral pelvic stabilisation along with enhanced trunk stability. All Jay products are designed to reduce fatigue and maximize function to provide clients the comfort and stability to go anywhere with confidence.

Here you find a selection of the Jay range of cushions and backs, which can be ordered via the Groove orderform.

For further Jay products, please refer to the Jay brochure.

JAY®

REHo20202

JAY® J² seat cushion

Ideal for clients who have moderate to more aggressive postural needs that may change over time, the Jay2 cushion is an ideal stable base of support.

The combination of the firm anatomically correct contoured base and fluid pads, provide excellent lateral stability as well as forward/rearward stability.



REHo20203

JAY® Combi seat cushion

Designed for clients who have a symmetrical posture needing minimum to moderate postural support.

The Jay® Combi features anatomically correct contours & deep leg troughs providing increased medial thigh support, whilst also providing trochanteric support to redistribute the load from the ischials



OPTIONS AND ACCESSORIES



JAY® J² back

Designed for clients who have minimal to moderate postural needs requiring taller trunk support, the Jay 2 back provides excellent posterior pelvic stabilisation and posterior thoracic support.

The J² back is also excellent at providing lateral support and addresses the need of clients with skin integrity issues.

REHo30003



JAY® J² deep contour back

Ideal for clients with moderate postural needs requiring more aggressive lateral support, the lightweight Jay deep contour back provides excellent pelvic and stabilisation and thoracic support, facilitating thoracic extension.

REHo30004

COMFORT SEATING

The Comfort Seating system provides clients with stylish seating that has been designed to maximise comfort whilst meeting mild to moderate requirements.

COMFORT SEATING

PROo2o100

Comfort seat pan

A firm padded board is a solid base for any seat cushion.

This seat pan with easy clean upholstery is available in the following sizes:

40cm, 44cm, 48cm wide
(15.5", 17", 19")

48cm & 55cm deep
(19", 21.5")



PROo2o101

Comfort seat pan with 5 cm (2") threshold

The comfort seat pan with added thresholds is designed for clients with their own cushion in mind.

The added thresholds give the benefit of stopping the movement between the cushion and the seat pan, giving a stable base for clients.



OPTIONS AND ACCESSORIES



Standard cushion - stretched rubber upholstery

This cushion has a 12 mm (1/2") beech wood base with a slightly contoured foam cushion, featuring a 2 cm (3/4") seat angle.

PROo2o2oo



Optimum comfort - stretched rubber upholstery

The Optimum comfort cushion includes a 12 mm (1/2") beech wood base with contoured foam cushion.

PROo2o2o1

COMFORT SEATING

PRO030000

Standard velour backrest - hard back

This anatomically shaped backrest features a luxury velour finish for added warmth and comfort and an easy clean rubber upholstery to the sides and rear.



PRO030001

Standard velour backrest - vertical elastic straps

This Standard Velour comfort backrest features flexible elastic straps instead of the hard wooden structure featured in the hard back version shown above.

The flexible elastic straps provide a softer and more comfortable ride whilst still providing adequate support.



Rear view



OPTIONS AND ACCESSORIES

PRO030002

Standard velour backrest - horizontal adjustable straps



This standard velour comfort backrest features tension adjustable straps which can be adjusted and shaped to the client's back, providing excellent support and flexibility.



Rear View

PRO030004

Deep contour backrest - passive



This comfort backrest features a deeper contour, anatomically shaped to provide support and stability across the entire back for clients with restricted movement.

Tension adjustable straps also provide additional support and flexibility for the client.



Rear view

PRO030003

Deep contour backrest - active



This comfort backrest features a deeper contour, anatomically shaped to provide support and stability around the lower back and is cut out around the shoulders so as not to restrict active clients.

Tension adjustable straps also provide additional support and flexibility for the client.



Rear view

PERFECT FIT SEATING SYSTEM

With versatile and comprehensive seating options, the Perfect Fit seating system is the ideal solution for any client requiring tailored Seating & Positioning.

Seating & Positioning adjustment is quick and easy with Perfect Fit. For the first time you can adjust the seat and back widths without any extra parts.

When it comes to the seating challenge, Perfect Fit won't take no for an answer.

PERFECT FIT SEATING

Perfect Fit seat & back

The Perfect Fit seating system is available in seat widths between 30cm - 61cm (12-24") and depths of 30cm - 56cm (12-22")

Seat and back width adjustment is possible in the following ranges;

Seat Width:	30cm - 43cm (12 - 17")	38cm - 51cm (15-20")	45cm - 61cm (19-24")
-------------	---------------------------	-------------------------	-------------------------

Back Width:	30cm - 33cm (12 -13")	36cm - 48cm (14-19")	43-61cm (17-24")
-------------	--------------------------	-------------------------	---------------------



PRFo20100

Perfect Fit contoured cushion 5 cm (2")

This 2" standard cushion is anatomically contoured to provide clients with added support and consists of a vented cover to improve airflow.

Available in all seat sizes.



PRFo30200

Perfect Fit Backrest Cushion

The Perfect Fit backrest cushion features 5cm (2") of adjustability built in along a vented and machine washable cover.

Available in all back sizes.



PRFo30400

Height adjustable push handles

The integrated push handles on the Perfect Fit seat system are detachable and height adjustable.



OPTIONS AND ACCESSORIES

POSITIONING ACCESSORIES



Lateral thoracic support - fixed

Available with 13cm x 13cm (5" x 5") or 15cm x 15cm (6" x 6") support pads

GRVo40300 -
GRVo40301



Lateral thoracic support - swing away

Swing your lateral support away with the click of a button, making side transfers easy.

Available with 13cm x 13cm (5" x 5") or 15cm x 15cm (6" x 6") support pads

GRVo40302-
GRVo40303



Wrap around lateral thoracic support - swing away

Similar to the swing away lateral thoracic, this support wraps around the front of the client providing improved stability and thoracic support

Available with 10cm x 10cm (4" x 4") 13cm x 13cm (5" x 5") or 15cm x 15cm (6" x 6") Support pads

GRVo40500-
GRVo40605



Swing away lateral knee support

Available with 10cm x 10cm (4" x 4") 13cm x 13cm (5" x 5") or 15cm x 15cm (6" x 6") support pads

GRVo40700 -
GRVo40805

RECARO® SEATING SYSTEM

Correct posture results in a pleasant journey. You'll arrive at your destination feeling relaxed and comfortable, but getting there in style is also important for many clients.

The range of Recaro seats available on the Quickie Groove tick all of the boxes. Stylish looks, adjustability and positioning and support where it really counts.

When it comes to the seating challenge, the Quickie Groove won't take no for an answer.

RECARO SEATING

RECo10002

Recaro® EXPERT M - manual recline

The Recaro Expert M provides clients with excellent seating and positioning with its extensive range of features.

The long and flat seat cushion makes it ideal for medium to larger clients. An adjustable lumbar support relieves stress on the spine. Seat extension provides better thigh support relieving muscle tension and the adjustable contoured back provides clients with excellent lateral support via adjustable side bolsters.



RECo10003

Recaro® EXPERT M - power recline

The Recaro Expert M - power recline provides the same features and benefits as the manual recline version with the added benefit of being able to recline the back via the hand control or most input devices.

Ideal for less able clients.



OPTIONS AND ACCESSORIES



Recaro® ERGOMED ES - power recline

No other seat offers you more comfort. The Ergomed ES is the top model of the Recaro product range.

The seat recline and adjustable side bolsters are both powered from the easily accessible control pad on the side of the seat, ensuring correct posture at all times.

RECo10004
RECo10005

Third party seat interface

If the seating systems we offer are not suitable, choose a third party seating interface.

This gives you the flexibility to add any seating system your client requires.

RECo10001



BACKREST RECLINE OPTIONS

The Quickie Groove has a recline option to suit everyone. Via tool adjustment, carer adjustment or client adjustment and via a manual or power solution.

REHAB & COMFORT SEATING

REHo40000

Manual adjustable backrest

The standard rehab manual adjustable back provides a fixed back angle from -5° to 12° .

Adjustment is still quick and easy. Simply remove the fixing bolts, set the desired angle and refit the bolts.



REHo40001

Manual recline - gas spring system

If your client requires a back that can be reclined whilst in the chair, then the gas spring system is for them.

With the client seated in their Groove, carers can recline this back from 0° to 30° with the press of a lever.



REHo40002

Power recline

Ideally suited for clients who want to adjust the back angle with no assistance from a Carer.

Powered from the hand control, the back can be reclined from -10° to 40° .



OPTIONS AND ACCESSORIES

PERFECT FIT SEATING



Manual recline backrest

Through Keipers, similar to those found on most cars, this option enables clients to recline their perfect fit backrest from -10° to 40° with the turn of a knob.

This assembly is hidden underneath the rear back shroud to give a slick and streamlined finish

PRFo40000



Power single pivot recline

Ideally suited for clients who want to adjust the back angle with no assistance from a Carer.

Powered from the hand control the back can be reclined from 0° to 40° .

This assembly is hidden underneath the rear back shroud to give a slick and streamlined finish

PRFo40001

SEAT TO BASE INTERFACE

The Groove power chair offers simple and versatile solutions for positioning.

A range of options are available. Manual preset tilt of 0 - 9 degrees, power tilt, power lift or a combination of both powered lift and tilt.

Simple modular units can be added or taken away as desired. This capability ensures that the Groove power chair is future-proof.

SEAT TO BASE INTERFACE

GRVo50000

Power tilt module

The Quickie Groove power tilt provides clients with 0° to 30° of tilt adjustment controlled through most input devices.



GRVo50001

Power lift module

The Quickie Groove power lift provides clients with 12" (30cm) of lift controlled through most input devices.



OPTIONS AND ACCESSORIES

POWER SEATING MODULE



Power lift & tilt combi module

The Quickie Groove lift & tilt module provides clients with 12" (30cm) of lift and 0° to 30° of tilt adjustment, powered through most input devices.

GRVo50002



Filler module

When a powered Seat Module is not required a filler module is used instead. This enables a 0°, 3°, 6°, or 9° seat angle to be set.

GRVo50003

ARMRESTS

The Quickie Groove offers a comprehensive range of armrests and accompanying sideguards to suit everyone. Once you have chosen your armrest we have a wide range of armpads from standard desk length pads to arm and hand pads designed specifically for Hemiplegic clients.

ARMRESTS

GRVo10001

Single post armrest & sideguard

(Standard Rehab, Comfort & Recaro Seating)

Available on Standard Rehab, Comfort and Recaro Seating, the detachable, height adjustable armrest can be fitted with a variety of pads and is available in short and long lengths.



GRVo10002

Flip-up armrest (Perfect Fit)

Only available on Perfect Fit seating the flip-up armrest makes preparing for side transfer quick and easy.

Choose from a variety of armrest pads



GRVo10003

Reclining armrest (Perfect Fit)

The Perfect Fit reclining armrest provides the ultimate positioning solution for clients who spend a large proportion of their time with the back reclined.

As the back reclines, the angle of the armrest adjusts automatically, reducing shear.

The armrest also detaches making side transfers possible.



OPTIONS AND ACCESSORIES

SKIRT GUARDS

Standard skirt guard

Available in short or long lengths this skirt guard is fixed to the armrest and can only be removed with the armrest.



Comfort skirt guard

For added comfort choose the Comfort Seating skirt guard with easy clean upholstery and soft padding.

Like the standard sideguard, this can only be removed with the armrest.



Perfect Fit skirt guard

The Perfect Fit skirtguard features removable, machine washable covers and a quick release mechanism, making detaching or moving the sideguard for side transfer quick and easy.



ARMPADS

GRVo10100 -
GRVo10101

Full length

36cm (14") long and padded for extra comfort, this is the tried and tested Quickie full length pad found on many Quickie chairs.



GRVo10102 -
GRVo10103

Desk length

25cm (10") long, ideal for clients who spend a lot of time behind a desk and padded for extra comfort.



GRVo10104 -
GRVo10105

Comfort armpads

The comfort armpad is contoured for extra stability and comfort.



GRVo10106 -
GRVo10111

Hemiplegic forearm pads

Available in three sizes the hemiplegic armpad cups the client's arm for extra stability and support.



OPTIONS AND ACCESSORIES

HEMI HANDPADS



Hemiplegic flat handpads

Available in two sizes and used in combination with the hemiplegic armpad, the wider profile gives clients extra hand support.

GRVo10200-
GRVo10203

COMFORT SEATING ARMPADS



Skai (vinyl) upholstery

The Comfort Seating armpad features an easy clean upholstery and extra padding.

PROo40000 -
PROo40003



P-Shaped skai upholstery

This Comfort Seating armpad has a wider profile.

PROo40004 -
PROo40007

LEGRESTS

Articulating, elevating and different angles of hangers, together with a wide range of footplates to suit every client preference.

CENTRE MOUNT LEGRESTS

GRVo20000-
GRVo20001

Fixed centre mount legrest

Ideally suited for clients requiring extra stability and support around their legs, the fixed centre mount legrest comes with calf pads and a flip up footplate and is available in 2 sizes:

Short: 20.4 cm - 31.6 cm (8" - 12.5")
Medium: 37 cm - 47 cm (13.5" - 18.5")
Long: 42 cm - 55 cm (16.5" - 18.5")



Powered centre mount legrest

Ideally suited for clients requiring elevation and support around their legs, the powered centre mount legrest is available in two sizes to suit different lower leg lengths.



Short: 25cm - 35cm (10" - 15") Long: 38cm - 48cm (14" - 19")

SWING-AWAY LEGREST RECEIVERS

Fixed width legrest receiver

The standard legrest receiver offering no width adjustment.



Width adjustable legrest receiver

This option is standard on Perfect Fit and Comfort Seating and available on Standard Rehab seating on 14", 16" and 18" seat frames.

The width adjustable legrest receivers now give clients the option of taking the hangers in or out from the seat



OPTIONS AND ACCESSORIES

SWING-AWAY LEGRESTS

Swing-away hangers - Push the button to swing away and detach. The range of swing away legrests is available in all frame colours.



90° hanger

GRVo20202



80° hanger

GRVo20201



70° hanger

GRVo20200

REHo20201



75° tapered hangers

The 75° tapered hanger with a 5cm (2") inset provides clients with additional support around their legs and makes manoeuvring through tight spaces easier than with the standard hangers.

GRVo20207

SWING-AWAY LEGRESTS

GRVo20203

Manual elevating legrest

The manual elevating legrest can be easily adjusted by the flip of a lever.

Swing away, detachable and complete with knee pads as standard.



GRVo20204

Powered elevating legrest

The powered elevating legrest is powered from the hand control and gives even the most challenged client the ability to adjust the legrest.

Swing-away, detachable and complete with knee pads as standard.



GRVo20205

Powered articulating legrest

Ideal for clients with severe clinical requirements and powered from the hand control, the articulating legrest automatically extends as it elevates.



OPTIONS AND ACCESSORIES

FOOTPLATES



Composite flip-up

Our tried and tested Quickie flip up footplate made from a composite plastic means it is lightweight and rigid.

GRVo30000



Aluminium angle adjustable

Made from aluminium this footplate will benefit clients looking for a more stable platform, and the ability to adjust the angle of the footplate.

GRVo30001



Aluminium flip-up

Only available on 90° hangers, the aluminium footplate can be easily flipped up for transfers and provides a stable and rigid support.

GRVo30002

FOOTPLATES (INC HEEL LOOP)

GRVo30003

Footboard

The one piece, angle adjustable footboard is constructed from lightweight aluminium giving clients the most stable platform available, whilst still retaining the benefit of detaching one side and flipping it up.



LEGREST ACCESSORIES

GRVo30100

Legrest strap



Heel loop

Standard on the composite flip up, angle adjustable and aluminium flip up footplates



OPTIONS AND ACCESSORIES

LEGREST ACCESSORIES



Standard calf pad

The padded calf pad cups client's legs whilst in the elevated position for added comfort and weight distribution.



Comfort seating calf pad

The comfort seating calf pad is made from a thick extra soft foam and an easy to clean cover for added comfort and convenience.

The Quickie Electronics control system breaks the mould with design, modularity and advanced performance, offering more functionality than any other control system on the market.

For the first time, dealers and therapists have the ability to programme and tailor drive and seating profiles to client needs, and assign multiple functions directly to hand control buttons making life easier for the user.

QC CONTROLS

DELo10000

QC-4 (Four button hand control)

Single profile hand control used for drive only options on the QC Control system.



DELo10001

QC-9 (Nine button hand control)

Single profile hand control used for drive, powered seating and lighting options on the QC control system.



OPTIONS AND ACCESSORIES

QR CONTROLS



QR-3 (Three button hand control)

Multiple profile hand control used for drive and powered seating options on the QR control system.

DELo20000



QR-7LT (Seven button hand control control lights and indicators)

Multiple profile hand control used for drive, powered seating and lighting options on the QR control system.

DELo20001



QR-2 (Two button compact joystick)

Compact multiple profile hand control used for drive and powered seating options. Best used in combination with a QR-ED display. Allows one fixed (power on / off) and one assignable optional switch input.

Also used for attendant and dual control.

DELo20003

DIRECT ACTUATOR CONTROL SWITCHES

DELo60001

Dual toggle switch

Two direction toggle switch allowing for remote operation of up to two powered seating functions. Allows seating to be operated simultaneously with driving.



DELo60002

Hex switch

Four direction switch joystick and two additional push button switches allowing remote operation of up to six powered seat functions. Allows seating to be operated simultaneously with driving.



OPTIONS AND ACCESSORIES



Two jack sockets for remote switches (heel)

Adaptable heel that can be retro-fitted to any QR hand control option, allowing two additional switch inputs, one fixed (power on / off) and one assignable optional switch input.

DELo60100



Speed potentiometer and toggle switch (heel)

Adaptable heel that can be retro-fitted to any QR hand control option, allowing remote on / off and assignable function via the toggle switch, and variable speed control on each drive profile using the potentiometer.

DELo60101



QR-ED Enhanced display module

Visual client display module that opens up the world of environmental control to the client. Allows customised client menus, access to Infra red audio / visual and environmental technology, as well as giving visual menu access on the QR system. Recommended with the QR-2 compact joysticks and all special control options.

DELo30000

CONTROL MOUNTINGS

In this section you will find all of the relevant control mountings for your Quickie Groove.

All control mounts shown here are compatible with every seating system. Please note a control mount receiver is required to fit each mount to your chair and will differ dependant on the options chosen.

QUICKIE HAND CONTROL MOUNTS

Fixed front mount

The fixed front mount is the Standard control mount for Quickie hand controls.



Swing-away front mount

Unlike the fixed front mount this Swing away mount offers clients horizontal adjustment of the hand control allowing users to swing it away and adjust it to their preference.



Fixed rear mount

Used to mount the hand control to the rear of your Quickie Groove for operation by carers or other third parties.



OPTIONS AND ACCESSORIES

Swing-away tray mounted centre control

The swing-away tray mounted centre control enables the hand control to be flipped over providing a flat surface when the tray is required and back around when the user requires access to the hand control.



Forus attendant control

The Forus attendant control enables clients to operate the Quickie Groove by pushing the adjustable handles in the desired direction providing more intuitive driving not achieved with the Standard control mount.



QUICKIE ENHANCED USER DISPLAY MOUNTS

Enhanced user display - fixed front mount

The fixed front mount is the standard control mount for the Quickie enhanced user display.

This mount provides angle adjustment through a series of easily adjustable knuckle joints.



Enhanced user display - swing away mount

The Swing away mount for the enhanced user display provides horizontal and angle adjustability, enabling clients to adjust their display.



Enhanced user display - flexi mount

The flexi mount is the most adjustable of all of the enhanced user display mounts, providing optimum positioning for viewing.



OPTIONS AND ACCESSORIES

Quickie mini joystick - fixed front & rear mount

This mount is used to mount the Quickie mini joystick to either the front or rear of the chair.

This mount provides angle adjustment through a series of knuckle joints.



QUICKIE MINI JOYSTICK MOUNTS

Quickie mini joystick - swing away front mount

Unlike the fixed front mount this swing away mount offers clients horizontal adjustment of the hand control allowing them to swing it away and adjust it to their preference.



CONTROL BOXES

QC motor controller

Quickie Consumer motor controller, capable of drive only operation.



QC-ACT-LT motor control

Quickie Consumer motor controller, operating up to two powered seating options plus lights and indicators.



QR motor controller

Quickie Rehab motor controller. Capable of drive, operating up to six powered seating options (when using QR-MAC), lights and indicators. Allows access to Intellidrive option on motors.



OPTIONS AND ACCESSORIES



Quickie multiple actuator control (QR-MAC)

Additional seating module on the QR system that upgrades the system from two to a maximum of six powered seating functions. Allows access to enhanced recline options.

DELo50000



Quickie environmental control (QR-ECM)

Additional "wired" environmental control module allowing access to nine additional switched inputs via the input device. Allows direct access to third party devices such as Dynavox communicators.

DELo50002



Quickie special control module (QR-SCM)

Allows access to all types of third party special control input devices. Operates switched and proportional devices, as well as scanning functions. Incorporates two stereo jack sockets, giving access to one fixed switch (power on / off) and three assignable switches.

DELo50001

SPECIAL CONTROLS

The Quickie Groove and Quickie electronics control system support a comprehensive range of Quickie and third party speciality controls to meet every client's need.

The Quickie electronics control system will also interface with any other speciality control system on the market.

SPECIAL CONTROLS

HMC mini joystick

HMC small movement joystick that connects via the QR-SCM module. Available as a hand control (DELo60201) or a chin control (DELo60203)



DELo60200

Five switch wafer board

Five switch input device that connects via the QR-SCM. Allows four directional switch driving plus a 5th optional mode switch.



DELo60202

Switch-It head array – adult

Third party switched head control that connects via the QR-SCM. Allows four directional switch driving (forward / reverse toggle) plus an optional mode switch.



DELo60203

Switch-It head array – paediatric

Third party switched head control that connects via the QR-SCM. Allows four directional switch driving (forward / reverse toggle) plus an optional mode switch.



OPTIONS AND ACCESSORIES



QR-2 chin control

Compact multiple profile chin control used for drive and powered seating options. Best used in combination with a QR-ED display. Allows one fixed (power on / off) and one assignable optional switch input.

DELo60300



Proportional chin control

Third party proportional chin control that connects via the QR-SCM. Allows 4 direction proportional driving.

DELo60301



Sip & puff control

Breath control option that connects via the QR-SCM. Allows driving in one of three different breath control configurations.

DELo60302



Switch-It proportional head control.

Third party proportional head control that connects via the QR-SCM. Allows four direction proportional driving (forward / reverse toggle) plus an optional mode switch.

DELo60204

The Power³

014867
015200

Quickie IQ QC/QR controls kit
Specialty Controls kit

2 Comprehensive kits of control parts to facilitate any client assessments or product simulation.



HEAD CONTROL MOUNTS

Powered swing away arm
- for chin and sip and puff controls



This option is used to mount all chin controls and can be operated by either a buddy button or ribbon switch, shown below.

DELo60311

Buddy button

Mounted on an adjustable fixed arm, users can control the powered swing away arm for the chin or sip and puff control by applying pressure to the buddy button.



OPTIONS AND ACCESSORIES



Ribbon Switch

Mounted on an adjustable fixed arm, users can control the powered swing away arm for the chin and sip and puff controls with the slightest touch to the ribbon switch.

DELo60312



Manual swing-away arm

This option is used to mount all chin controls. You can also mount buddy buttons directly onto the arm to make operating the Groove easier.

DELo60310



Buddy button switches....

Red
Green
Yellow
Blue

DELo60400
DELo60402
DELo60404
DELo60406

WHEELS

The Quickie Groove offers an extensive range of drive wheels and castors to suit most types of drive terrain and environment. Choose from a variety of tyre treads and sizes, available in solid, pneumatic or OKO filled (puncture resistant).

DRIVE WHEELS

GRVo80000
GRVo80001

13" Pneumatic
13" Solid

This slick tread is ideal for indoor and outdoor use available in solid and pneumatic.



GRVo80002
GRVo80005

14" V groove pneumatic
14" V groove pneumatic
c/w OKO fluid

Like the 13" drive wheels the V groove is ideal for indoor and outdoor, available in pneumatic or OKO filled for puncture resistance.



GRVo80003
GRVo80004

14" Knobby-pneumatic
14" Knobby-solid

The Knobby tyre really comes into its element when exposed to the elements! Ideally suited for outdoor use giving extra grip in rough terrain.

Available in Solid or Pneumatic on 10 or 13kph motors.



OPTIONS AND ACCESSORIES

CASTOR WHEELS



8" x 2" Solid

GRVo8o101



9" x 2.5" Solid
9" x 2.5" Pneumatic

GRVo8o102
GRVo8o103



10" x 3" Solid
10" x 3" Pneumatic

GRVo8o104
GRVo8o105



Suspension castor system

Tackling rough terrain in the countryside or for city living, the castor suspension system will handle it all.

Configured to the client's weight the castor suspension systems provide a softer, more comfortable ride.

GRVo8o201 -
GRVo8o203

BATTERIES & CHARGERS

The Quickie Groove offers a choice of batteries and chargers.

Clients who are likely to use the Groove for drive only and basic functions will find the 60Ah batteries suitable, but if the Groove is equipped with features like an enhance user display, lift or tilt and lights and indicators then the larger 72Ah battery would be more suitable.

BATTERIES

GRVo90100

60 Ah maintenance free



GRVo90101

72 Ah maintenance free



OPTIONS AND ACCESSORIES



220V / 8A automatic charger UK

GRVo90200

220V / 8A automatic charger
Continental

GRVo90201

For use with 60Ah batteries



220V / 10A automatic charger UK

GRVo90203

220V / 10A automatic charger
Continental

GRVo90204

For use with 72Ah batteries and 13kph
motors.

ACCESSORIES

The Quickie Groove offers an extensive range of accessories that have been designed to assist users with thier day to day living.

GRVo90500

Crutch holder



GRVo90501

Rear view mirror



GRVo90505

Serial PS/2 infra red computer mouse port receiver

Plugged directly into any PC or laptop this accessory allows clients to control their computer's mouse using their Groove hand control.

(Enhanced user display also required)



OPTIONS AND ACCESSORIES



Curved headrest 20cm x 11cm
(8" x 4.5")

GRVo40000



Curved headrest 23cm x 11cm
(9" x 4.5")

GRVo40001



Contoured headrest

With adjustable wings the contoured headrest can be adjusted to give optimum support to clients

GRVo40002



Comfort Seating - standard headrest

This comfort seating headrest provides a soft velour finish for warmth and easy-clean rubber upholstery.

GRVo40003



Comfort Seating - flexible headrest

This Comfort Seating headrest provides the same benefits as the standard headrest with the added ability to adjust the flexible wings around the client's head for added support

GRVo40004

GRVo40100

5 cm (2") positioning strap



GRVo40101

5 cm (2") padded positioning strap



GRVo90504

24V to 12V dc/dc converter

This accessory converts the 24V output to 12V enabling clients to run 12V equipment such as mobile phone chargers.




GRVo90504

Kerb climber

Climbing kerbs up to 10cm (4") high is quick and easy with the Groove kerb climber.

Only available on Groove RWD.

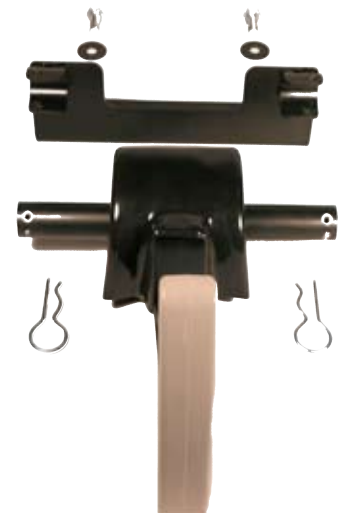

OPTIONS AND ACCESSORIES

The Quickie Groove offers an extensive range of joystick knobs for clients with restricted activity or who are just looking to personalise their Hand Control.

The functionality and size of each knob has been designed for finger, palm, chin or hand operation, reducing the effort needed to operate the Groove.

JOYSTICK KNOBS

Golf ball joystick knob



DELog0001

T-bar joystick knob



DELog0002

Soft ball joystick knob



DELog0003

Ball joystick knob



DELog0004

Chin cup joystick knob



DELog0005

Stick joystick knob



DELog0006

DELog0007

Mushroom joystick knob



Swing-away tray table

The Quickie Groove tray table provides clients an easy to clean and stable surface to work or eat from. Once finished simply swing out of the way or detach and store away.

Available in the version shown below or with the hand control mounted in the centre of the tray (refer to parts manual).



OPTIONS AND ACCESSORIES

For the first time Quickie dealers are able to easily reconfigure and personalise the Quickie Groove using the revolutionary Quickie IQ software kit and handheld programmer.

QUICKIE IQ ACCESSORIES



Hand held programmer (HHP)

The Quickie hand held programmer enables the dealer to adjust basic settings and profiles on a client's Groove.

The Quickie hand held Programmer is available exclusively to approved Quickie Dealers

GRV090502



Quickie IQ programming software kit.

The Quickie IQ Software Kit enables the dealer to reconfigure and programme a client's Quickie Groove allowing drive parameters and more advanced settings to be adjusted.

Also required for firmware / software upgrades and loading / saving files.

GRV090503



Web IQ is the world's first interactive online programming tool, enabling dealers to create, review and store unique customer profiles in a safe and secure environment.

Features include:

- Ability to download chair configuration and reconfigure the Quickie Groove
- Review chair history
- Save drive profiles on a secure server (located at Sunrise Medical)
- Download files to a chair using Quickie IQ.

www.sunmedwebiq.com


LIGHTS

Choose from our standard STVZO approved lights and indicators or exclusive to the Quickie Groove, LED lights and indicators. These project a more focused and intense beam providing greater visibility. The LED Lights also draw less power from the battery and have a considerably longer life then the standard lights.

STANDARD LIGHTS

GRVo60000

Standard front & rear lights and indicators
(StVZO Approved)

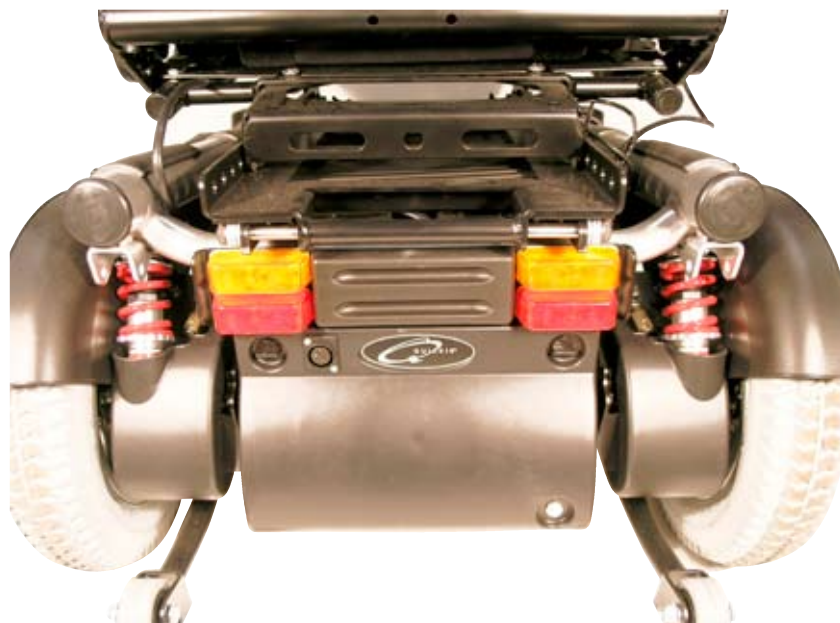

OPTIONS AND ACCESSORIES

LED LIGHTS

LED front & rear lights and indicators

GRVo60001





Groove™

QUICKIE
iq
ELECTRONICS



Sunrise Medical Limited
High Street Wollaston
West Midlands
DY8 4PS, UK

STLPS™
Sunrise Training and Education Programmes

KNOWLEDGE FOR THE FUTURE



Tel: +44 (0) 1384 44 66 88
www.quickiepower.com
For international enquiries e-mail:
International@sunrisemedical.de

SR0052

**INNOVATE.
BUILD FOR THE FUTURE**